

MYH7 antibody



Catalog No: #AT24727

Host Species: Mouse

Clonality: Monoclonal antibody (mAb)

Species Reactivity:
Human; Mouse; Rat;
Rabbit; Avian

SDS-PAGE MW:
223 kDa

Applications:
Western blotting
Immunofluorescence
Immunohistochemistry

Package Size:
#AT24727S 20 μ l
#AT24727M 50 μ l

Description

Product Name MYH7 antibody

Purification Monoclonal antibody is purified by affinity isolation.

Specificity MYH7 monoclonal antibody detects endogenous level of total MYH7 protein.

Immunogen Type
Recombinant Protein

Immunogen Description
Recombinant Protein

Target Name MYH7

Other Names CMH1; MPD1; SPMD; SPMM; CMD1S; MYHCB

Concentration 1.0 mg/ml

Storage Store at -20°C for long term preservation (recommended). Store at 4°C for short term use.

Note In order to get the best results using various techniques and preparations, we recommend detecting optimal working dilutions by titration.


Application Details

Western blotting: 1:250~1:1000 (recommended 1:500)

Immunohistochemistry (IHC): 1:50~1:200

Immunofluorescence (IF): 1:50~1:200

Images

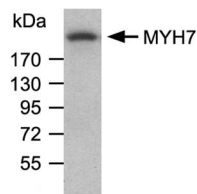

Western blot analysis of extracts from human skeletal muscle tissues, using MYH7 mouse monoclonal antibody #AT24727.



Tel: 025-85390309; 85300038

Orders: sales@sciben.com

Support: support@sciben.com

Website: www.sciben.com